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ELECTRIC HEATING PAD**INTENDED USE**

The electric heating pad with a maximum power of 60 W can be used in situations where it is necessary to warm a certain part of the body. It has a soothing power and relaxes the muscles, which helps with arthritis, neuralgia and rheumatism. The temperature controller allows you to adjust the temperature in 3 steps. The heating pad is for home use only and is not suitable for use in hospitals, doctors' offices, etc.

Characteristics:

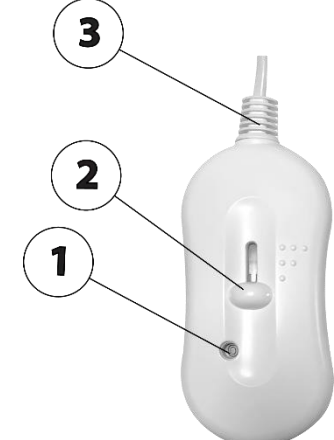
- Soothing power of heat for joints, nerves and muscles
- It is recommended for treatment of arthritis and rheumatic diseases
- Soothes neuralgia, relaxes muscles and improves blood circulation
- Protection against overheating
- Temperature control using the remote control on the power cord
- Fast and even heat distribution

Use the device only for the prescribed purposes. Any other use is considered a case of misuse. The user/operator and not the manufacturer will be responsible for any damage or injury caused by this misuse. Please note that this device is not designed for commercial or industrial use. The warranty will not be valid if the device is used for commercial, industrial or similar purposes.

TECHNICAL PARAMETERS

	VOLTAGE / FREQUENCY	230 V / 50 Hz
	POWER INPUT	60 W
	COLOR	White
	MATERIAL	100% Polyester
	POWER CABLE LENGTH	1.5 + 0.5 m
	DIMENSIONS	150 x 80 cm
	WEIGHT	1 kg

PRODUCT PARTS

			
1	Heating pad	1	LED indicator
2	Operated with power cord	2	Switch for levels
		3	Power cable

EXPLANATION OF SYMBOLS	
	Read the instructions for use.
	Do not dispose of in normal household waste. Instead, turn to recycling centers in an environmentally acceptable way. Please take care to protect the environment.
	The product complies with the applicable European directives and a conformity assessment method of these directives has been carried out.
	Protection by double insulation II. No grounding is required.
	The heating pad is not intended for children under 3 years of age. For children over 3 years of age, the heating pad may only be used under adult supervision.
	Do not use the heating pad if it is bent or folded.
	Do not insert any sharp objects, such as safety pins, into the heating pad; they could damage the electrical cables.
	Do not dry the heating pad in the dryer.
	Do not clean the heating pad by any professional cleaning method.
	The heating pad can be washed at 30° without damaging it. Gentle washing without spinning. It is recommended to wring the heating pad only by hand and, it is necessary to reduce the amount of laundry in the washing machine.
	Bleaching the heating pad is prohibited.
	It is forbidden to iron the heating pad.

SAFETY INSTRUCTIONS FOR HEATING PADS

- **Read all safety warnings and instructions. Failure to follow safety warnings and instructions may cause electric shock, fire and/or severe injury.** Keep all safety warnings and instructions carefully for future reference.
- Carefully unpack the heating pad, taking care not to throw away any part of the packaging material before you have found all the components of the product.
- The heating pad may be used by children over 8 years of age, as well as persons with reduced physical, sensory or mental abilities or lack of experience and knowledge, provided that they are supervised or have been instructed in the safe use of the heating pad and understand the hazards involved. Otherwise, there is a risk of faulty operation and injury. Children must not play with the heating pad. Children must not carry out cleaning or user maintenance without supervision.
- Never use a heating pad if your hands or feet are wet.
- Use and store the heating pad away from flammable and volatile substances. Never place the heating pad near heaters, open flames, or other appliances or equipment that are sources of heat.
- Do not expose the heating pad to extreme temperatures, direct sunlight, excessive humidity or place it in excessively dusty environments.
- The heating pad is designed for use at home. Do not use it in industrial environments or outdoors!
- The heating pad must not be used as a pet cushion. There is a danger of biting electrocution, with possible death of the animal as a result.
- Before connecting the heating pad to an electrical outlet, make sure that the voltage indicated on the product label matches the voltage at the electrical outlet.
- Do not route the power cord near hot surfaces or over sharp objects! Place it so that it is not accidentally pulled out of the socket, and so that no-one can trip over it!

- Do not stick any sharp objects, such as safety pins, into the heating pad; they could damage the electrical cables. Do not sew the heating pad with a needle!
- The heating pad must not be rolled or covered during operation. Never fold the heating pad when it is switched on. Always make sure that, if used on a reclining bed, the heating pad is laid out flat and cannot be rolled up or the cord caught or pinched. If it cannot be guaranteed that the above situations can occur, do not use the heating pad on a reclining bed.
- Do not cover the heating pad or use it as a bedsheet.
- Do not place the heating pad under the mattress.
- Make sure that the heating pad remote control is not covered by blankets, cushions or similar objects and is within easy reach.
- Prolonged use of the heating pad at a high level can cause burns on the skin. The heating pad should not be used by disabled persons, children or persons sensitive to heat.
- **WARNING!** The user should not fall asleep when using the heating pad.
- When the heating pad is spread out on the bed, do not let children jump on it. If they do, check thoroughly that every part of the pad is in order.
- Never disassemble or repair the heating pad yourself. Tampering with the equipment will void the warranty.
- If the electric heating pad is no longer to be used, it must be turned off and disconnected from the power source by pulling the plug of the power cable from the power outlet. Never pull the heating pad out of the socket by the power cable while it is switched on and if your hands are wet or damp. Make sure the power cord is not damaged.
- In case of any damage to the power cable, immediately switch off the heating pad and disconnect it from the power source. The power cord cannot be replaced. The heating pad must be disposed of in accordance with the applicable regulations.
- Store the heating pad folded (or twisted) not too tightly, without bending the heating element in a dry place. Do not place heavy objects on the pad or twist it. After storage, carefully inspect the pad before further use. Do not use the pad with a damaged power cord.
- If dirt appears on the pad, follow the instructions for cleaning the heating pad in the "Cleaning and Maintenance" section.

OPERATION AND ASSEMBLY

INSTRUCTIONS FOR USE

- Make sure the heating pad fits your mattress and then lay it flat directly on the bed mattress to leave room for a pillow. The edge where there is no power cord with the control belongs to the legs. The heating pad should be firmly attached to the mattress. Make sure that the heating pad is not wrinkled or kinked.
- Make the bed.
- The pad is designed so that it can warm even the uncovered part of the bed.
- The pad should not be tucked into the sides of the bed.
- Connect the power cord with the remote controller to the heating pad, and then connect the power cord plug to an electrical outlet (make sure the power cord is not damaged or tangled in any way).
- When the heating pad is in use, the remote controller must not be covered. Make sure it is not placed under the heating pad, cushion or other objects, or near other heat sources, it could overheat and be damaged.
- Push the switch to position I, II or III. The LED indicator will light up, and the heating pad will start to heat up.
- To stop heating, push the switch to position 0.

OVERHEATING PROTECTION

- The heating pad is equipped with a safety system that protects the pad from overheating and damage.
- In the event of a malfunction, the system switches the pad off and cannot be used again. In this case, an authorised service center must be contacted.

(*) Authorized service: the manufacturer's or importer's service department or other person who is qualified,

authorized and competent to carry out this type of repair. The heating pad must be handed over to such a service in the event of repair.

CLEANING AND MAINTENANCE

- If the heating pad becomes dirty, switch off the appliance and pull the power cord plug out of the power outlet.
- Disconnect the power cord with the controller from the heating pad.
- Small dirt can be removed with a damp cloth in lukewarm water.
- Remove larger scale dirt by washing. Place the heating pad in the washing machine and wash gently in lukewarm water, using a wool detergent. Do not put any sharp objects that could damage the heating pad into the washing machine with the heating pad.
- To ensure maximum durability, we recommend washing the heating pad a maximum of 3 times.
- After cleaning the dirty area, do not use a hairdryer to dry the heating pad, but leave it loosely draped to dry. The heating pad should not be allowed to dry draped.
- Do not expose the heating pad to direct sunlight.
- Never turn the heating pad on to dry it.

ENVIRONMENTAL PROTECTION

The crossed-out wheeled-bin symbol on the products or accompanying documents means that used electrical and electronic products should not be mixed with general household waste. For proper disposal, recovery and recycling, please deliver these products to designated collection points, where they will be accepted free of charge. Alternatively, in some countries, you may return your products to your local retailer when purchasing an equivalent new product. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in penalties in accordance with national regulations.

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EN: This product is guaranteed for 24 months from date of sale or removal from storage. All manufacturing defects found during the warranty period will be repaired at no charge. To submit a request for warranty repairs, fill out the warranty card legibly and completely and attach it to the device. Customers pay the cost of transport when sending a device for repairs. Keep the product's original packaging in a safe place.

The warranty does not cover:

- Damage caused to the device during transit or from improper storage
- Malfunctions caused by incorrect operation or maintenance
- Malfunctions caused by product and material wear
- Malfunctions caused by using the device for purposes other than those intended
- Damage caused by tampering or modifications
- Any incomplete work not discovered on the product when it was sold

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